





BASIC BIOLOGICAL MICROSCOPE FOR SMALL LABS



Motic's microscope series Panthera S introduces a new quality and performance level to basic routine and scientific work. This microscope combines distinct improvements in image quality with new illumination features. The Panthera S is equipped with the redesigned SC Plan Achromatic Infinity optics. These objectives combined with N-WF widefield eyepieces with a 20mm Field of View deliver un-

matched image quality and definition. The Panthera S stand features a reversed 4-fold nosepiece with a traditional x/y stage with a large travel range, and the illumination is provided by an LED light source in a Fixed Koehler setup. Enhanced objective performance, long-lasting illumination, and a user-friendly design - the Panthera S is perfectly tailored for modern education and small lab environments.



360° rotatable head (storage saving space)



Quadruple nosepiece Mechanical stage

Integrated USB for camera power & USB for battery pack


EFFICIENT LED ILLUMINATION SETUP

The Panthera S is equipped with a low consumption LED light source in a Fixed Koehler setup without field diaphragm. The light intensity indicator displays the light intensity as a sequence of LEDs around the revolving nosepiece. With a USB port in the stand, Panthera S can connect to a power bank and function independently from the main power supply. The standard Abbe condenser with NA 0.9/1.25 with aperture diaphragm and slider slot is ready to carry sliders for Phase contrast and Darkfield.





Motic's Panthera S features a rackless stage with a standard slide holder, as well as an upper limit stop, eliminating the risk of breaking the glass slide while focusing. The minimum increment of 2 microns in z-movement facilitates smooth focusing of samples under high magnifications. The standard eyepiece tubes can be adjusted to accommodate different physiques. In addition, all tubes are rotatable, a helpful space saving feature for storage.

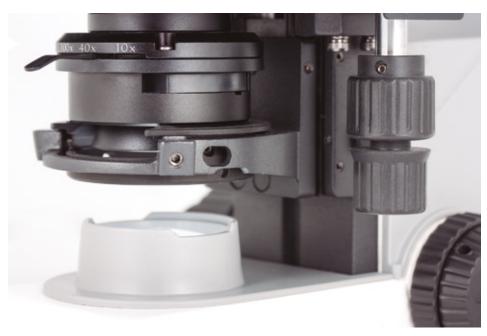

Our moticam range offers the choice between several imaging options. Besides the standard USB cameras, we offer solutions like the Moticam Full HD cameras, tablet solutions, and Wi-Fi cameras working with our free "MotiConnect" App. Know more at www.moticamseries.com

NEW SC OBJECTIVES AND EYEPIECES

For Panthera S, Motic redesigned the CCIS® Infinity optics to create a series of new Plan SC Achromatic objectives. These lenses present a flat image for a 20mm Field of View with optimal color fidelity and contrast. The new 100X/0.8 "Dry objective" is especially favorable in environments where the time-consuming immersion procedure on a "normal" 100X objective is considered too difficult and complicated.





SAMPLE CHARACTERISTICS

If desired, Panthera S can be upgraded with accessories for different contrast methods. A Phase contrast slider solution for unstained samples from botany, zoology, and medicine, a Darkfield slider, and a POL contrast option for birefringent samples.

PANTHERA S



Optical system	Colour Corrected Infinity Optical System (CCIS®)
Observation system	Binocular/Trinocular Siedentopf type 25° inclined, 360° rotatable (storage saving), 360° swivelling
Interpupillary distance (mm)	48-75mm
Trinocular light split	50/50 fixed
Eyepieces	N-WF10X/20
Eyepieces diopter adjustment	+/- 5 dpt
Nosepiece	Reversed quadruple
Objectives Plan SC	4X/0.1 (WD 15.5mm), 10X/0.25 (WD 17.4mm), 40X/0.65/S (WD 0.6mm), 100X/1.25/S-0il (WD 0.14mm)
Objective mounting thread	W 4/5" x 1/36" (RMS standard)
Stage	Mechanical stage
Stage size (mm)	140x135
Travel range (mm)	75x50
Upper limit stop	Preset; adjustable
Condenser	N.A. 0.90 / 1.25 Abbe condenser with slider slot
Focus mechanism	Coaxial; tension adjustment
Minimum fine focus precision (μm)	2
Z-axis movement (mm)	25
Illumination type	0.5W LED, Fixed Koehler
Illumination features	Nosepiece LED light intensity indicator
Connectivity	USB for external camera power & USB for battery pack
Power supply	110-240V (CE)
Dimensions (mm)	282x210x395 (Binocular), 282x210x396 (Trinocular)
Contrast techniques	
Brightfield	Yes
Phase contrast	Slider
Simple Polarization	Yes
Simple Darkfield	Slider
LED Fluorescence	Yes - Intermediate module

DISCOVER THE COMPLETE PANTHERA SERIES AT www.moticpanthera.com

www.moticeurope.com

WETZLAR, GERMANY

SALES OFFICE

Christian-Kremp-Straße 11, 35578 Wetzlar, Germany T. +49 6441 21001 0 | Fax: 49-6441-210 0122

BARCELONA, SPAIN

SALES OFFICE & EUROPEAN DISTRIBUTION CENTRE

C. Les Corts 12, Pol. Ind. Les Corts 08349 Cabrera de Mar, Barcelona, Spain T. +34 93 756 62 86 l Fax: 34 93 756 62 87







