

Moticam BTU10

TABLET SOLUTIONS

The Moticam BTU10 is a tablet microscopy solution, consisting of a customized 10" Android tablet connected to a microscopy imaging camera.

The android tablet is easily detachable and replaceable and is completely separate from the camera, therefore providing maximum flexibility.

The Moticam BTU10 will help microscope users visualize, edit and share the microscope images, enabling computer-free digital microscopy.

www.moticamseries.com | www.moticeurope.com

10" TABLET MICROSCOPY SOLUTION





Moticam **BTU10** TECHNICAL SPECIFICATIONS

Camera Specifications

Sensor Type	CMOS
Sensor Size	1/2,5"
Live resolution	1280 x 720 pixels, 1280 x 960 & 800 x 600
Image capture resolution	5.0MP
Video capture resolution	720P (HD)
Data transfer	USB
Pixel Size	2.2µm x 2.2µm
Lens Mount	C-Mount
Tablet bracket holder	For tablets from 7" to 10"

Tablet Specifications

Size	263 x 173 x 9.5 mm
System	Android 5.1
Pre-installed App	MotiConnect
Language	Multi-Language
CPU	1.9GHz Quad core
RAM	2GB DDR3
NAND FLASH	16GB
Display	10.1" (IPS) 16:10 Ratio
Ports	MicroSD/MMC/TF Card (Support up 32 GB), 3.5mm Audio, DC 5V Power Input, Mini HDMI output
Camera	Front (2M Pixels)
Support Display Format	MPEG-1/2 (DAT, MPG, VOB, TS), MPEG-4, H.263 (3GP, MOV, MP4), H.264 (AVI, MKV, MP4, MOV, 3GP), RM, RMVB, FLV, WEBM
Support Audio Format	MP3, WMA, MP2, OGG, AAC, M4A, OGA, FLAC, APE, MP1, WAV
Picture Format	JPG, JPEG, BMP, PNG
Package Includes	MicroSD 16Gb included, HDMI Cable, USB 2.0 Cable, Motic Images Plus 3.0, Motic 4-dot calibration slide, Charger 100V-240V (CE)


MOTIC IMAGES PLUS 3.0 SOFTWARE & MOTICONNECT APP

THE [AllinOneBOX] CONCEPT

Everything that you may need to work with your Moticam is included in one box

Motic Images Plus 3.0

By connecting your Moticam to your computer you will be able to work with our well-known software, that comes included in the package. View, capture, edit, measure, make reports... all its standard features and the new ones packed in a new user-friendly interface.

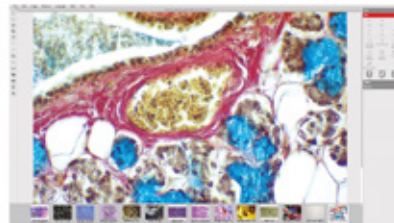

MotiConnect App

View, edit, capture and store your images in the blink of an eye, view live images on multiple wireless devices and share your work with anyone in just a few clicks. The MotiConnect App makes data sharing possible in a collaborative environment, encouraging teamwork and collective learning experiences.





Canada | USA | Europe | China



www.moticmicroscopes.com

Motic Instruments, Inc. (Canada)

130-4611 Viking Way, Richmond, BC V6V 2K9 Canada
Tel: 1-877-977-4717 | Fax: 1-604-303 9043
info@motic-america.com

Motic Instruments USA Inc.

6508 Tri-County Parkway Schertz, TX 78154 USA
Tel: 1-800-275-3716

Motic Europe (Spain)

C. Les Corts 12, Pol. Ind. Les Corts. 08349 Cabrera de Mar, Barcelona, Spain
Tel: 34-93-756 6286 | Fax: 34-93-756 6287

Motic Incorporation Ltd. (Hong Kong)

Unit 2002, L20, Tower 2, Enterprise Sq. 5, 38 Wang Chiu Rd, Kowloon Bay, Kowloon
Tel: 852-2837 0888 | Fax: 852-2882 2792

Design Change: The manufacturer reserves the right to make changes in instrument design in accordance with scientific and mechanical progress, without notice and without obligation.

Motic Incorporation Limited Copyright © 2002-2019. All Rights Reserved.

Updated: October 2019

