

**#153054**  
w/Sage #153048


**#153054B** w/Powder #153044



White  
**#153010**



Green  
**#153015**

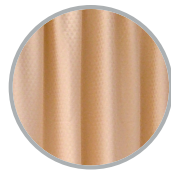

Apricot  
**#153020**




Oatmeal  
**#153028**



Marine  
**#153029**



Norway  
**#153035**

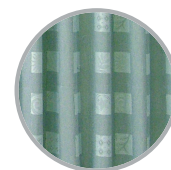


Sage  
**#153048**



Rouge  
**#153040**



Autumn  
**#153042**

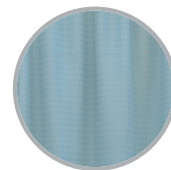

Powder  
**#153044**



Frost  
**#153045**



Lavender  
**#153046**

**Screen Frames and Panels are Sold separately**

### Beamatic™ Screen Frame #153054 - #153054B

The 4 sections of the screen are attached to each other by a riveted hinge which allows the frame to smoothly fold and unfold as needed. The S-shaped base of each screen section on the frame also increases stability for use in busy environments. Privacy screen panels, either in economical vinyl or attractive designer fabrics are needed to complete the screen. One panel is required for each section and is secured by top and bottom snap in rods. Panels are discounted when purchased with the screen.

#### Features:

- Available in 2, 3, 4 or 5 section models
- Four sections measures 68”H x 73”W
- Requires one vinyl or fabric panel per section
- Anodized aluminum construction provides a lightweight and durable solution eliminating the possibility of paint chipping and rust
- 5 year warranty on Frames
- 90 Days on Screen Panels